

# Inion FixOn™

Biocomposite Suture Anchor

**INION**



# Inion FixOn™ Biocomposite Suture Anchor



The fully threaded Inion FixOn™ Biocomposite Suture Anchor provides excellent performance over the healing period in rotator cuff surgery.

Inion FixOn™ Biocomposite Suture Anchor is made of state-of-the-art biocomposite material containing PLGA and osteoconductive  $\beta$ -TCP. Based on *in vitro* degradation testing, the polymer component of the device is expected to fully degrade within 2-3 years. The products are provided sterile and preloaded on a practical, single-use driver with two or three nonabsorbable sutures.



## Benefits

- Excellent visibility of the implant during the operation**  
 Inion FixOn™ Biocomposite Suture Anchor is the only colored biocomposite suture anchor on the market. The pigment (D&C Green No. 6) has routinely been used in implants and sutures for decades.
- Highly visible, strong sutures (white/striped/blue)**  
 Suprabraid™ sutures are prepared from ultra-high molecular weight polyethylene (USP #2, metric 5).
- Fast insertion and improved fixation strength over the healing period**  
 The cannulated, fully threaded design with an additional thread in the proximal end increases cortical fixation.
- Awl-only insertion technique**  
 The high-strength implant material and optimized instrumentation enable the awl-only technique when preparing the insertion site in most bone densities.
- Optimal material composition**  
 Biocompatibility, strength and degradation properties are supported by test data and clinical studies of similar material composition.

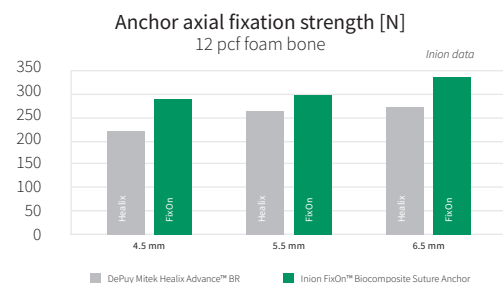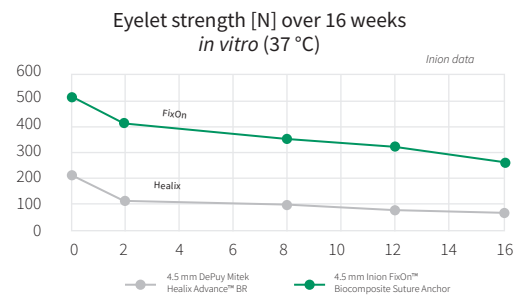


## Ordering Information

### Implants:


- FOB-4502** 4.5 mm Suture Anchor with two sutures  
**FOB-4503** 4.5 mm Suture Anchor with three sutures  
**FOB-5502** 5.5 mm Suture Anchor with two sutures  
**FOB-5503** 5.5 mm Suture Anchor with three sutures  
**FOB-6502** 6.5 mm Suture Anchor with two sutures  
**FOB-6503** 6.5 mm Suture Anchor with three sutures

### Instruments:

- INS-9630** Reusable tap for 4.5 mm Suture Anchors  
**INS-9631** Reusable tap for 5.5 mm Suture Anchors  
**INS-9632** Reusable tap for 6.5 mm Suture Anchors  
**INS-9633** Reusable awl for 4.5 mm Suture Anchors  
**INS-9634** Reusable awl for 5.5/6.5 mm Suture Anchors

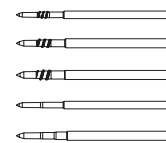

- INS-9640** 2nd level tray  
 (INS-9640 can be accommodated in INS-9619 Inion® Shoulder instrument tray)

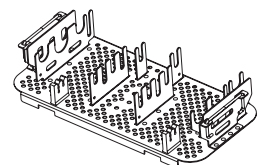

**INION OY**  
 Lääkäriinkatu 2  
 FI-33520 Tampere, FINLAND  
 tel: +358-10-830 6600  
 email: info@inion.com  
 internet: www.inion.com

**INION INC**  
 2800 Glades Circle  
 Suite 138, Weston  
 FL 33327  
 USA  
 Toll-free tel: 866-INION-US  
 tel: 954-659-9224-fax: 954-659-7997